

Gness 150 CS 📼

A RELIABLE PARTNER

UNMATCHED SIMPLICITY IN ROCKWELL HARDNESS TESTING





THE PERFECT INTRODUCTION TO ROCKWELL HARDNESS TESTING

AT HOME EVERYWHERE: FOR INCOMING GOODS, WORKSHOPS, LABORATORIES AND SERIAL PRODUCTION

LED TEST SPACE LIGHTING Built-in LED workspace lighting simplifies the precise positioning of the test piece, thus enhancing operating convenience.

AFFORDABLE PRICING

SOPHISTICATED CONSTRUCTION, ULTRA-SIMPLE OPERATION



COMPACT DESIGN -LATEST TECHNOLOGY

Modern features provided as standard:

- I Controlled force application (1 kg to 250 kg)
- I Direct depth-gauge system (0.05 μm resolution)
- Robust welded steel frame and sheet-steel covers
- I 7" color touchscreen operation with Qpix TE software
- I Digital data delivery via USB interface



TEST TABLE HEIGHT ADJUSTMENT

Via stable, ultra-precise roller-bearing spindle guide. Solid, nomaintenance structure. All devices are available with a Ø 25 mm table mount (optional $\frac{3}{4}$ " adapter available) facilitating the use of a wide range of test tables and devices.



MADE IN AUSTRIA

Unmatched in its class! We know how important a long-serving superior-quality device can be and guarantee excellent quality from our Qness 150 CS ECO series. Developed and manufactured in Austria!

SUPPORTED TEST METHODS

\bigcirc	ROCKWELL DIN EN ISO 6508, ASTM E-18							
	HRA	HRB	HRC	HRD	HRE	HRF		
	HRG	HRH	HRK	HRL	HRM	HRP		
	HRR	HRS	HRV	HR 15-N/T/W/X/Y				
	HR 30-N/T/W/X/Y			HR 45-N/T/W/X/Y				

BRINELL

HBT (not acc. to standards)




PLASTICS TESTING

49.03 N 132.9 N 357.9 N 961 N

Ų

 \rightarrow DIN EN ISO 18265, DIN EN ISO 50150, ASTM E140

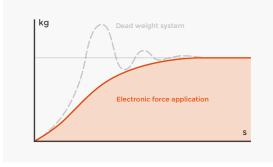

FULLY AUTOMATED TEST CYCLE

Electronic force application and closed-loop control

ESS

O C K W E L L

NORM-COMPLIANT WITH A DIRECT DEPTH-GAUGING SYSTEM

2020 10:24:52

91.01.2020 10.32 46

31.01.2020 11.02.38

28,79

HINC

Start test

QATM configures the Qness 150 CS ECO with a direct, ultra-precise HEIDENHAIN depth-gauging system, positioned exactly on the axis of the indenter. This prevents measurement inaccuracies caused by deformation in the system.

Another unbeatable benefit: Indenting depth can be gauged directly and can be easily calibrated in accordance with the latest Rockwell norms.

SOPHISTICATED OPERATION

INTUITIVE FEATURES AND OPERATION



START BUTTON

Start button at the front of the tester ensures a simple and rapid start to hardness testing (particularly when wearing gloves), and there's also touchscreen control.



OPERATION DISPLAY AND QPIX TE SOFTWARE

QATM has developed a 7" color display with a robust sheet steel frame for the Qness 150 CS ECO. The display uses capacitivetouch technology and has been optimized to guarantee operational precision and rapid reaction times.

The new Qpix TE software guarantees simple operation with larger buttons and a more modern and intuitive software interface.

Features:

- I Extra-large, clearly-visible presentation of test results
- I High-/low-quality evaluation of results featured in green and red according to the tolerance levels set
- I Test results list for 999 values
- I Statistics overview (min/max/range/cp/cpk/average)
- Conversions
- I Surface corrections
- I Simple test method changeover
- I Wide range of operating languages
- I Password protected area for expanded range of settings



REMOVABLE DOWNHOLDER

No need for long retooling breaks due to inaccessible test points. The clamping head can be removed via 2 screws. Now, even test points in hard-to-reach contours can be accessed with ease.



DIGITAL DATA EXPORT

Test results and all related information can be transferred to other data storage mediums via the USB port at the side.

TESTING

HARDNESS

OCKWELL

READY-TO-USE

COMPLETE PACKAGE

QNESS 150 CS ECO HARDNESS TESTING DEVICE

DAkkS & ASTM-conformant Rockwell indenter included
Hardened test anvil

ADVANTAGES

Rapid deliveryFactory-calibrated

INDIVIDUAL OPTIONS

- I Test blocks
- Additional indenters
- I Machine pedestal
- I Clamping devices from the comprehensive QATM accessories program
- I Customer-specific clamping mechanisms



+



Height adjustment	manual / spindle			
Test height	250 mm			
Throat depth	157 mm			
Test table	Ø 100 mm			
Max. workpiece weight	100 kg			
Weight of basic device	77 kg			
Test sequence	fully automated / electronic force application			
Software	Qpix TE			
Display	7" capacitive touch (color)			
Ports	1x USB			
Options	Additional indenters, hardness test blocks, machine pedestal, prisms			

1 - 250 kg (9.81 - 2452 N)



QNESS 150 CS ECO PEDESTAL

No suitable table in the testing area? No problem! The Qness 150 CS ECO matches perfectly with the machine's optional pedestal. The superior-quality subframe is color-matched to the hardness testing device. The machine pedestal includes a spacious, lockable compartment in which to store accessories. Beechwood (40 mm thick) with a spindle bore.

Dimensions:

B 750 x T 700 x H 840 mm

Max. surface load 1250 kg

ONLINE PRODUCT-CONFIGURATOR

Test force range

For more equipment and accessories go to the online product configurator at www.qatm.com



Online konfigurátor >



07



ATM	Qness	GmbH	

Emil-Reinert-Str. 2 57636 Mammelzen Germany Reitbauernweg 26 5440 Golling Austria

Phone: +43 6244 34393

ATM Qness GmbH

Phone: +49 2681 9539 0 Fax: +49 2681 9539 27

PREMIUM QUALITY

Fax: +43 6244 34393 30

info@qatm.com www.qatm.com





Verder Scientific is a business field belonging to the Verder Group and sets standards in the development, manufacture and sale of laboratory and analytics devices. Used in quality control, research and development for test-piece preparation and the analysis of solids.

For several decades our companies have supplied production plants and research institutes, laboratories for quality testing and analytics, all kinds of technical specialists and scientists with modern, reliable devices to solve the many and varied challenges they face. CARBÓLITE

ELTRA Retsch

VERDER SCIENTIFIC

SCIENCE FOR SOLIDS